

AD-A040 788

ARMY INST OF DENTAL RESEARCH WASHINGTON D C
MODIFIED FILM HOLDER FOR ENDODONTICS.(U)
JUN 77 O E PAQUETTE, R O SEGALL, C E DEL RIO

F/G 6/12

UNCLASSIFIED

1 OF 1
ADA
040788



NL

END

DATE
FILMED
7-77

ADA 040788

REPORT DOCUMENTATION PAGE		READ INSTRUCTIONS BEFORE COMPLETING FORM
1. REPORT NUMBER (6)	2. GOVT ACCESSION NO. (12)	3. RECIPIENT'S CATALOG NUMBER
4. TITLE (and Subtitle) Modified Film Holder For Endodontics		5. TYPE OF REPORT & PERIOD COVERED Submission of Paper. Jan 77 - June 77
7. AUTHOR(s) (10) Omer E./Paquette, Ronald O./Segall Carlos E./del RIO		6. PERFORMING ORG. REPORT NUMBER
9. PERFORMING ORGANIZATION NAME AND ADDRESS US Army Institute of Dental Research Walter Reed Army Medical Center Washington, DC 20012		8. CONTRACT OR GRANT NUMBER(s)
11. CONTROLLING OFFICE NAME AND ADDRESS US Army Medical Research and Development Command ATTN: (SGRD-RP) Washington, DC 20012		10. PROGRAM ELEMENT, PROJECT, TASK AREA & WORK UNIT NUMBERS (16) Program Element 62110A Project No. 3A162110A825 Task Area No. 00-Work Unit No.
14. MONITORING AGENCY NAME & ADDRESS (If different from Controlling Office) (9) Rept. For Jan - Jun 77 (12) 9p.		12. REPORT DATE (12) 10 June 77 117
		13. NUMBER OF PAGES 6
		15. SECURITY CLASS. (of this report) UNCLASSIFIED
		15a. DECLASSIFICATION/DOWNGRADING SCHEDULE
16. DISTRIBUTION STATEMENT (of this Report) This document has been approved for public release and sale; its distribution is unlimited.		
17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report) DDC RECEIVED JUN 22 1977 A		
18. SUPPLEMENTARY NOTES		
19. KEY WORDS (Continue on reverse side if necessary and identify by block number) SNAP-A-RAY II, film holder, modification.		
20. ABSTRACT (Continue on reverse side if necessary and identify by block number) Proper film positioning and stabilization in spite of projecting instruments, rubber dam clamps, and frame have always been a challenge, particularly for students and assistants using the long cone, paralleling technique. Simple modifications to the new, Rinn Snap-A-Ray II film holder and positioner virtually eliminate the usual difficulties.		

DDC FILE COPY

DD FORM 1 JAN 73 1473 EDITION OF 1 NOV 65 IS OBSOLETE

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

038670

LB

ACCESSION NO.	
WHD	Write Section <input checked="" type="checkbox"/>
UDG	Unit Section <input type="checkbox"/>
UNANNOUNCED	<input type="checkbox"/>
JUSTIFICATION	
BY	
DISCUSSION (AVAILABILITY CODE)	
DATE	WHD, WUD, or UG
A	

MODIFIED FILM HOLDER FOR ENDODONTICS

OMER E. PAQUETTE, A.B., D.D.S.
 RONALD O. SEGALL, D.M.D., M.S.
 CARLOS E. del RIO, B.S., D.D.S.

- LTC, DC Division of Clinical Research, US Army Institute of
 Dental Research, Walter Reed Army Medical Center,
 Washington, DC 20012
- MAJ, DC Endodontic Resident, US Army Institute of Dental
 Research, Walter Reed Army Medical Center, Washington,
 DC 20012
- COL, DC Chief, Division of Clinical Research; Director,
 Endodontic Residency, US Army Institute of Dental
 Research, Walter Reed Army Medical Center, Washington,
 DC 20012

Properly angulated, well centered radiographs during bio-mechanical instrumentation are important for consistent quality results in endodontic therapy. Proper film positioning and stabilization in spite of projecting instruments, rubber dam clamps, and frame have always been a challenge, particularly for students and assistants using the long cone, paralleling technique. Simple modifications to the new, Rinn Snap-A-Ray II film holder* and positioner virtually eliminate the usual difficulties.

The radiolucent plastic, Snap-A-Ray II has an anterior, slip-in film holder at one end, and a posterior, hemostat-like film holder at the other, and is designed for use in taking bisecting, periapical exposures using either the short or long cone (Figure 1A). The posterior film holder, has a broad, thick jaw on which the patient bites to stabilize the film. The film holder positions and stabilizes the film packet, while a stainless steel rod connects the film holder with the ring-like aiming device which is placed flush with the end of the cone when the exposure is made. As supplied, this rod has a fixed, 15⁰ bend, standardizing the angle to be bisected. Two types of rod are provided; one with a right angle bend in addition to the 15⁰ bend, for use with posterior exposures, and the other, a straight rod having only the 15⁰ bend, for use in the anterior (Figure 2A).

*Rinn Corporation, 1212 Abbott Drive, Elgin, Illinois 60120

The broad jaws of the posterior film holder prevented correct film positioning when used with the rubber dam, and so they were modified. Excess plastic was ground away, the surfaces smoothed and then polished with pumice on the office lathe (Figure 1). To modify Snap-A-Ray for use with the paralleling technique, the 15⁰ bend in both the anterior and posterior rods were removed (Figure 2).



Using the modified Snap-A-Ray (Figures 3 & 4), the patient stabilizes the film packet parallel to the long axis of the tooth by holding the free end of the film holder. Mesially angled posterior exposures are easily managed by appropriately adjusting the film position in the holder distally. With this adjustment the long cone may be angled so that the rod does not interfere. With the modified Snap-A-Ray and a radiolucent, rubber dam frame, excellent radiographs are routinely obtained with no need to disturb placement of the rubber dam, a boon to experienced and student practitioner, alike.

REQUEST FOR REPRINTS:

Doctor Omer E. Paquette
Division of Clinical Research
US Army Institute of Dental Research
Walter Reed Army Medical Center
Washington, DC 20012

CAPTIONS

- FIGURE 1. The Snap-A-Ray II film holder, (A) before and (B) after modification.
- FIGURE 2. The anterior and posterior metal rods (A) before and (B) after modification. The upper rod in both groups is the anterior.
- FIGURE 3. The Snap-A-Ray II fully assembled for posterior use, (A) before and (B) after modification. Note the straightened rod at (B) which has had the 15° bend seen in (A) removed.
- FIGURE 4. The Snap-A-Ray fully assembled for anterior use, (A) before and (B) after modification.

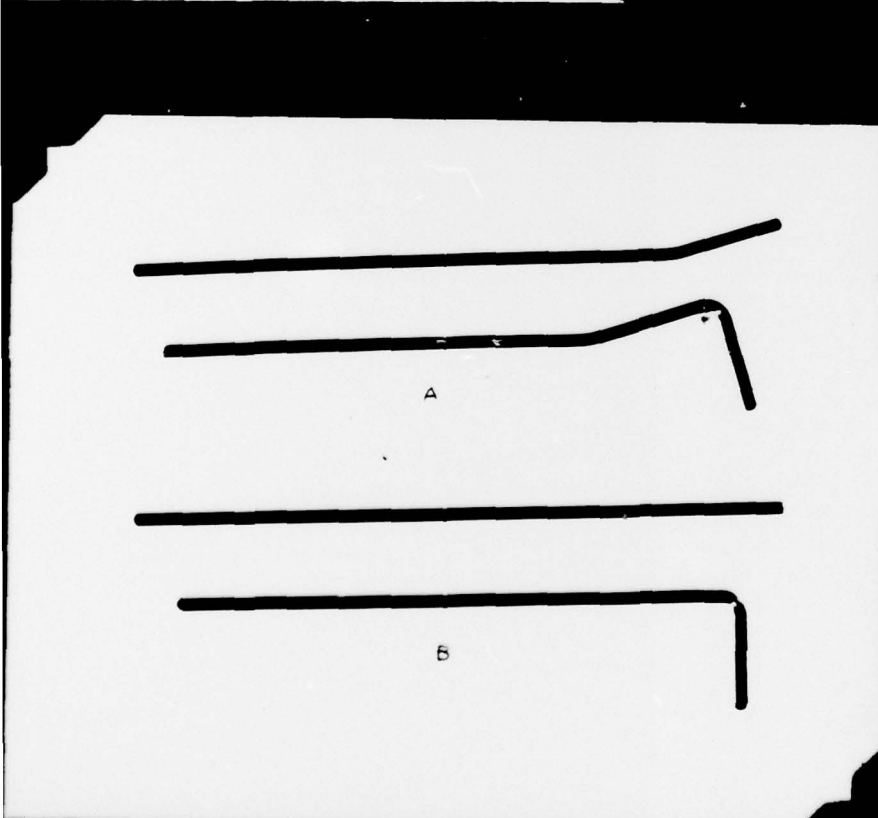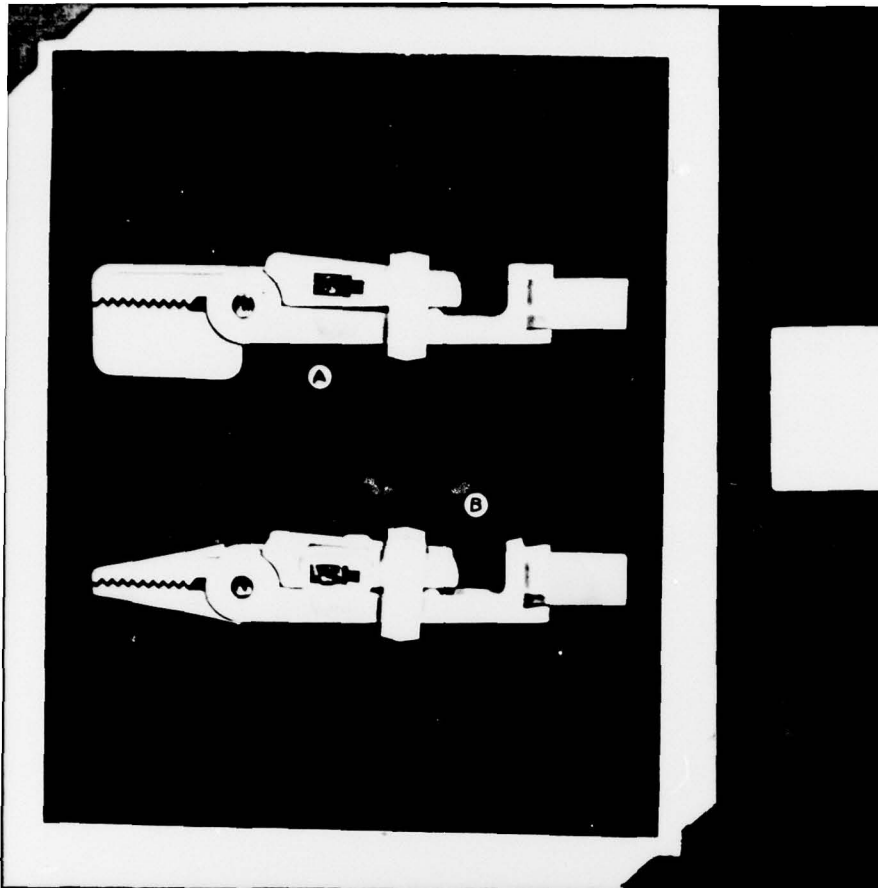

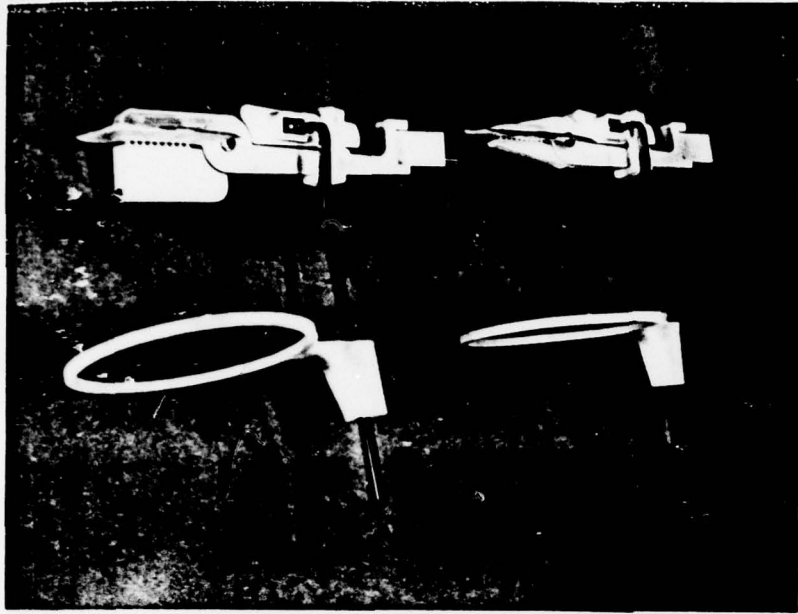


FIGURE 3

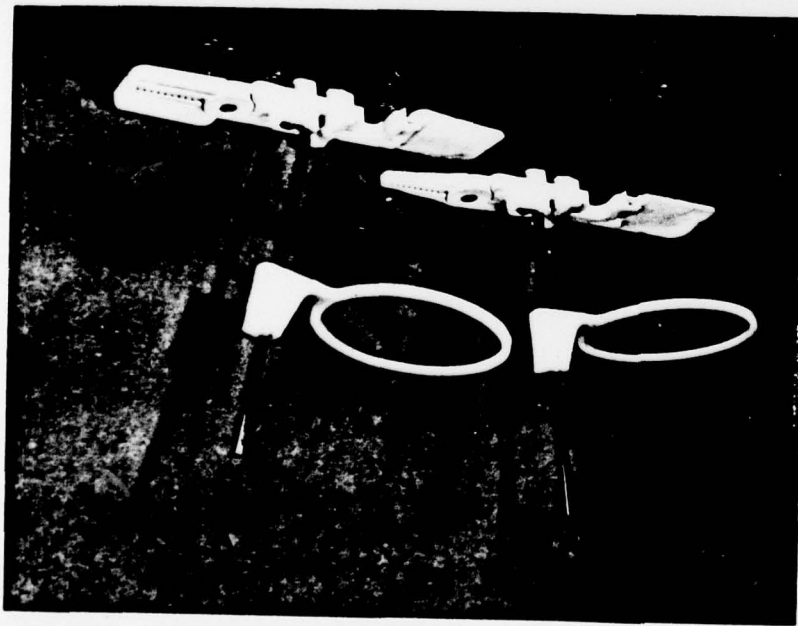


FIGURE 4

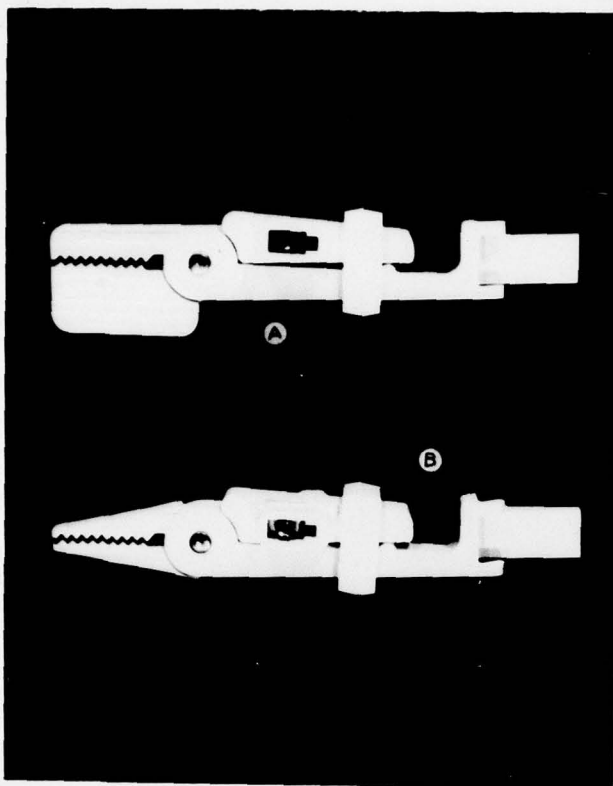


FIGURE 1

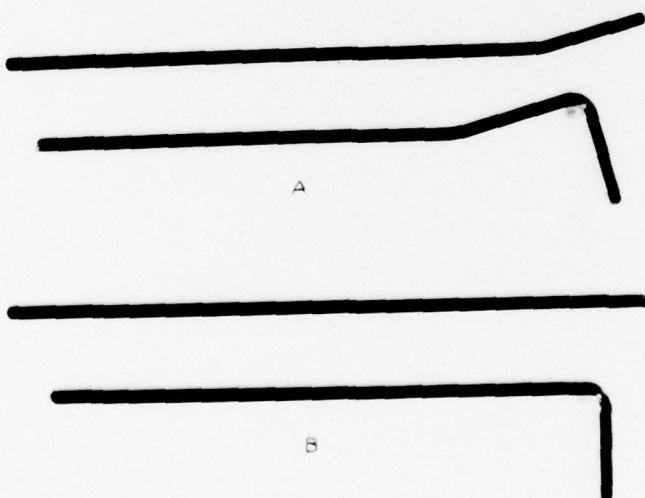


FIGURE 2

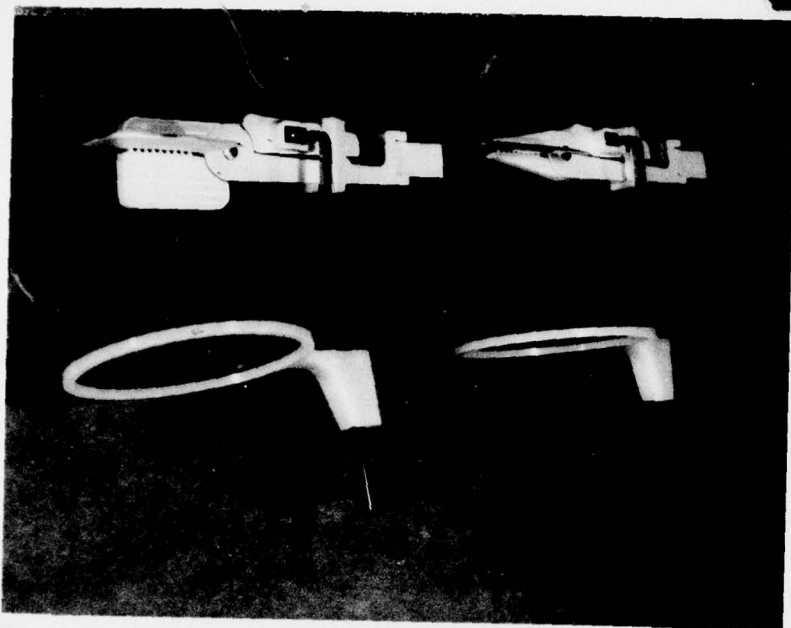


FIGURE 3

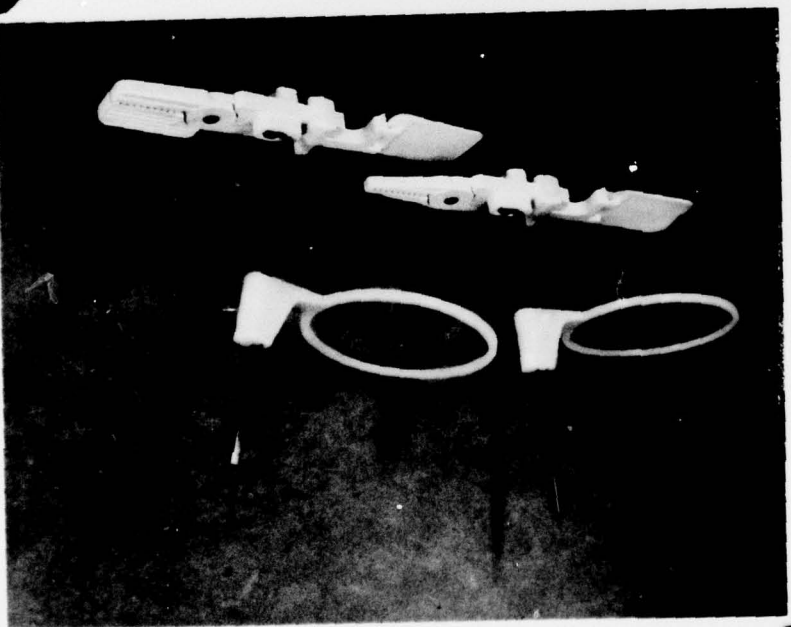


FIGURE 4